

SPECIFICATIONS There are 17 standard frequency options for the following products: 72.1, 72.2, 72.3, 72.4, 72.5, 72.6, 72.7, 72.8, 72.9, 74.7, 75.3, 75.4, 75.5, 75.6, 75.7, 75.8, 75.9 MHz

TRANSMITTERS

| ST-300 17-CHANNEL TH | RANSMITTER (CATALOG NO. 71912-000) |
|----------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| RF Frequency Range | 72 to 76 MHz |
| Modulation: | FM ±25 KHz deviation |
| Signal-To-Noise Ratio | 58 dB (77 dB with E.D.R.) |
| Maximum Deviation | ±25 KHz |
| Audio Input | Balanced XLR-3F and unbalanced 1/4" |
| Antenna | 1/4" wave telescoping omnidirectional |
| Audio Controls | Audio input level, monitor jack volume, E.D.R. |
| Power Requirements | 13-20 Vdc or 12 Vac; 115Vac 60 Hz @ 300 mA plug-in wallpack power supply |
| Dimensions | H 1-3/4" (4.5 cm) x W 7-1/2" (19.2 cm) x D 5.9" (15 cm) |
| FCC ID | B5DM522 |
| Backlit LCD Display | Ch., Freq, Audio Meter, Input setup, RF level, and E.D.R. |
| RF Output Power @ 3M Normal High | 25K μV/m 80K μV/m |
| PST-170 17-CHANNEL T | RANSMITTER (CATALOG NO. 71914-000) |
| Audio Input Characteristics | Mic TA4M connector, 7.75 mV RMS input for ±25 KHz deviation. Impedance, 10K Ohms, nominal. Aux. Input, Female 3.5mm. 100mV Unbalanced, 100K Ohm nominal |
| Antenna | Permanently attached trailing wire |
| Modulation: | FM ±25 KHz deviation |
| Frequency Response (System) | 100 Hz to 10,000 Hz |

58 dB (77 dB with E.D.R.)

(2) AA batteries, Alkaline, NiMH or NICAD

8 hours Alkaline, 8 hours NiMH, 4 hours NICAD

115 µ seconds

B5DM521

In accordance with Federal Communications Commission (FCC) rules and regulations, Telex SoundMate personal listening systems are intended only for use in providing auditory assistance to the handicapped and are not intended for use by the general public.


80K µV/m @ 3 m

SR-400 17-CHANNEL RECEIVER (CATALOG NO. 71913-000) (2) AA batteries, Alkaline, NiMH or NICAD 20 hours Alkaline, 16 hours NiMH, 8 hours NICAD **Power Requirements** Audio Frequency Response <3 dB Variation (100 Hz-10 KHz) Modulation FM ±25 KHz deviation Sensitivity 0.5 µV typical, 1.0 µV maximum, 12 dB SINAD Image Rejection >65 dB Signal-To-Noise Ratio >77 dB with E.D.R. <2% T.H.D. Distortion into 8 Ohms: 2.0V / 15mW 3.0V / 80mV into 32 Ohms: 2.0V / 10mW 3.0V / 50mW Audio Output @ 10% distortion Antenna Type omnidirectional in earphone cord Volume, On/Off, High frequency boost control, E.D.R. Audio Controls Earphone Connector Type 3.5 mm stereo or mono H 4" (101 mm) x W 2-3/4" (70 mm) x D 1" (25 mm) Dimensions Backlit LCD Channel, Hi Freq, E.D.R. **Visual Indicators** Bottom of unit Charger Contacts FCC Complies with Part 15 SR-50 RECEIVER (CATALOG NO. 71155-XXX*) (2) AA batteries, Alkaline, NiMH or NICAD 30 hours Alkaline, 20 hours NiMH, 10 hours NICAD **Power Requirements** Audio Frequency Response <3 dB Variation (100 Hz-10 KHz) Sensitivity 0.5 µV typical, 1.0 µV maximum, 12 dB SINAD Signal-To-Noise @ 1 mV Input >60 dB <2% T.H.D. Distortion into 8 Ohms: 2.0V / 15mW 3.0V / 80mV into 32 Ohms: 2.0V / 10mW 3.0V / 50mW Audio Output @ 10% distortion into 8 Ohms: omnidirectional in earphone cord Antenna Type Controls Volume/On/Off Switch Dimensions H 4" (101 mm) x W 2-3/4" (70 mm) x D 1" (25 mm) Visual Indicators Red LED On Indicator Earphone Connector Type 3.5 mm stereo or mono Charger Contacts Bottom of unit FCC Complies with Part 15

RECEIVERS

* The last three digits in the transmitter catalog number (XXX) will be determined by the frequency that is selected.



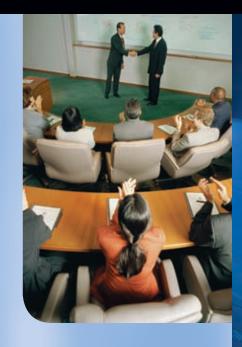

TELEX



System Signal-To-Noise

Pre-emphasis

FCC ID

Maximum RF Power

Power Requirements

 Americas
 Telex Communications Inc. 12000 Portland Ave South, Burnsville, MN 55337, USA

 USA
 Phone: 1-800-392-3497, Fax: 1-800-955-6831

 Canada
 Phone: 1-866-505-5551, Fax: 1-866-336-8467

 Latin America
 Phone: 1-952-887-5532, Fax: 1-952-736-4212

©2004 TELEX COMMUNICATIONS, INC. SPECIFICATIONS SUBJECT TO CHANGE WITHOUT NOTICE. PA-20652 USA 6/04



listening systems







Personal Listening System

Telex SoundMate[™] personal listening systems help overcome background noise and poor building acoustics that can make listening difficult for the hearing impaired listener. A Telex base station, portable belt pack transmitter, your choice of receivers, and a wide assortment of accessories allow you to meet the needs of individuals who require hearing assistance. These listening systems are even compatible with a variety of hearing instrument styles for hearing aid users. Let Telex, the acknowledged leader in personal assistive listening products, help you select the system that's right for any listening environment.



AUDIO SOURCE: Amplifier, Direct Microphone, or Mixer.

System Versatility - A sound system amplifier, mixer or wired microphone can plug directly into the ST-300 base station transmitter. The ST-300 then sends the audio signal to an unlimited number of personal receivers.

TRANSMITTERS



FRONT

ST 300

ST-300 TRANSMITTER The ST-300 Base Transmitter features 17 user selectable frequencies digitally controlled from the front panel control buttons, a headphone jack with adjustable level for input signal monitoring, and a backlit LCD display that shows Channel, Frequency and Audio Input Level. On the back of the unit is a balanced XLR-3F as well as an unbalanced 1/4" input. Front panel controls and LCD display allow the user to select mic, line, and 70 volt input options, attenuate the input levels and change hi/low RF transmit power. A new feature called Enhanced Dynamic Range (E.D.R.) works with the SR-400 receiver to provide clearer, more dynamic audio.





THE PST-170 BELT PACK

TRANSMITTER The PST-170 Belt Pack Transmitter may be set to any one of the 17 channels which is displayed on the backlit LCD, and operates on frequencies in the 72-76 MHz band. It is lightweight, battery operated and includes a small lapel microphone for convenient portable use. The PST-170 has two audio input jacks. The Auxiliary input is designed to allow the use of audio devices such as teacher's aids or tape players. The second jack is normally used with a microphone and is located on the top panel. The PST-170 is also compatible with the E.D.R. feature of the SR-400.



THE SR-50 SINGLE CHANNEL RECEIVER provides clear

reception and the simplicity of a fixed channel. This economical receiver operates on one of 17 fixed narrow-band frequencies in the 72-76 MHz band. Ergonomic raised volume control knob make level adjustments easily accessible by feel. Recessed headphone jack provides extra protection for earphone connections. Two AA batteries give up to 30 hours continuous battery life.

SYSTEMS



SM-2 SYSTEM Includes:

(1) ST-300 base transmitter, (4) SR-50 receivers (4) SEB-1 earbuds (1) SoundMate wall plaque. (SM case sold separately).

SMP-2 SYSTEM Includes: (1) PST-170 portable belt pack transmitter (6) SR-50 single channel receivers (6) HED-2 collapsible headphones (14) AA Alkaline batteries (1) Deluxe system carrying case (1) HM-2 head worn microphone (1) User manual for SR-50 and PST-170

R E C E I V E R S

THE SR-400 17 CHANNEL RECEIVER is the perfect solution to multiple transmitter hearing assistance systems. The SR-400 features an advanced digital PLL synthesizer to tune all 17 narrow-band frequencies offered in the 72-76 MHz band and a backlit LCD display. The receiver also includes a special high frequency contour filter and boost switch to increase the intelligibility of the audio signal, and E.D.R. for increased dynamic range. Additional features include a battery saving automatic power shutoff when the earphone is unplugged.





FRONT

BATTERY CHARGER

THE BC-100 BATTERY CHARGER

charges two SR-50's, SR-400's or PST-170's or any combination of transmitters and receivers.